Orgonzapper Standard, Orgonzapper El Silverado, Clarkzapper

Operating instructions and safety information

Application:

The zapper is designed as a wellness device for relieving muscle tension. It has not been tested or approved for the treatment of illnesses.

The zapper can be worn anywhere on the arms or legs using the elastic strap. It can also be placed on any other part of the body without the strap. The electrodes should be in light contact with the skin.

How it works:

The zapper generates a weak 'tingling current' of 9V pulsed direct current. This is perceived by many users as relaxing and invigorating.

Safety:

The zapper operates with extremely low currents of maximum 20mA, which means that there is no risk of danger. Nevertheless, it should not be used by people with pacemakers or during pregnancy.

Inserting and replacing the battery:

For flight safety reasons, the zapper is delivered without a battery. A 9V PP3 battery is required. To insert the battery, loosen the cover with your fingernail at the corresponding slots and lift it off. The battery must lie flat on the bottom of the device without any cables underneath, otherwise the cover will not close. The push button must also be positioned on the side facing away from the battery. Ensure that the battery is inserted with the correct polarity.

Switching on and off:

The red LED light flashes when the zapper is switched on and the battery is installed correctly.

Care and maintenance:

Discolouration of the skin electrodes is normal. Oxidation occurs, particularly with the silver electrodes, which you will recognise as 'tarnishing' from silver cutlery. As this is pulsed direct current, this will always be more pronounced on one of the two electrodes. There is no need to constantly polish the electrodes to a shine. It is best to wipe the electrodes with a damp cloth and diluted lemon juice. Metal polishes are not necessarily skin-friendly and should only be used in extreme cases. Then remove all residues with a damp cloth.

Please be gentle with the small slide switch, and the zapper will serve you well for many years.

Warranty:

The manufacturer guarantees the trouble-free functioning of the device for 6 months from receipt, provided it is used properly.

Disposal:

Due to its electronic components, the zapper must be disposed of as electronic waste when it reaches the end of its service life. Used batteries should be disposed of at appropriate collection points.

Manufacturer

Orgone Power (Pty) Ltd
166 Frances Street
Observatory 2198
Gauteng
Südafrika
info@orgoniseafrica.com
Tel +27102370180

Importer for the EU and authorised representative

Orgone Power Europe UG (haftunbgsbeschraenkt)
Dudenstrasse 57
10965 Berlin
Deutschland
kat@orgoniseafrica.com
Tel +4917670712371

EG - Declaration of conformity

The manufacturer

Orgone Power (pty) Ltd
166 Frances Street
Observatory 2198
South Africa
info@orgoniseafrica.com
+27 10 237 0180

hereby declares that the products

Orgonzapper Standard
Orgonzapper El Silverado
Essential Clarkzapper
Zapper bracelet

comply with the provisions of the following directives:

General Product Safety 2001/95/EU
Electromagnetic Compatibility (EMC) 2004/108/EC
Low Voltage Directive 2014/35/EU
ROHS 2011/65/EU


The articles comply with the following standards required for the CE mark:

EN 55014-1:2006+A1:2009+A2:2011
EN 55014-2:1997+A1:2001+A2:2008

Authorised representative for technical documentation in accordance with Directive 2006/42/EC:

Orgone Power Europe UG (limited liability)
Dudenstrasse 57
10965 Berlin
ph: +49 176 707 12371
email: kat@orgoniseafrica.com

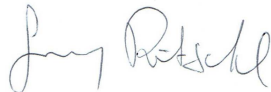Place and date of issue: Observatory, Johannesburg, 1 February 2024

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Orgone Power (pty) LTD